

The OrthoSave® Comp provides controlled iceless compression therapy for the reduction of swelling and pain. In addition, the OrthoSave® DVT offers mobile, mechanical blood clot prevention. With certain conditions where a physician prescribes heat or contrast therapy, OrthoSave® Comp also delivers great outcomes.

YOUR DAILY PRESCRIPTION *(default prescription in parentheses below)*

Use for \_\_\_\_\_ (30) minutes with 1 hour off between sessions \_\_\_\_\_ (3-5) times per day.


1 WHAT'S IN THE BOX?

- ☒ OrthoSave® Comp
- ☒ OrthoSave® DVT
- ☒ Power cord
- ☒ Umbilical hose
- ☒ Cold & compression wrap
- ☒ DVT wraps *(if applicable)*
- ☒ Distilled water
- ☒ UPS return label

\*Note: The machine color may be beige or dark grey

2 INITIAL SET-UP

ORTHOSAVE® COMP

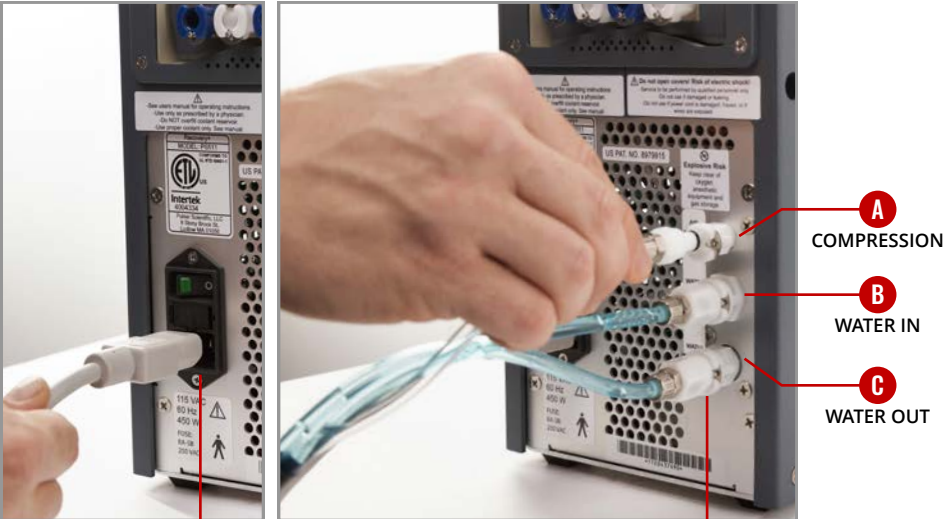

- STEP 1

Plug in power cord
- STEP 2

Connect umbilical hose lines to **A** COMPRESSION, **B** WATER IN and **C** WATER OUT

ORTHOSAVE® DVT



- STEP 3

Connect a *single DVT wrap* to **D** CALF LEFT


COMPLETED SET UP



3 STARTING THERAPY

- STEP 1

Apply all wraps prescribed by your physician.

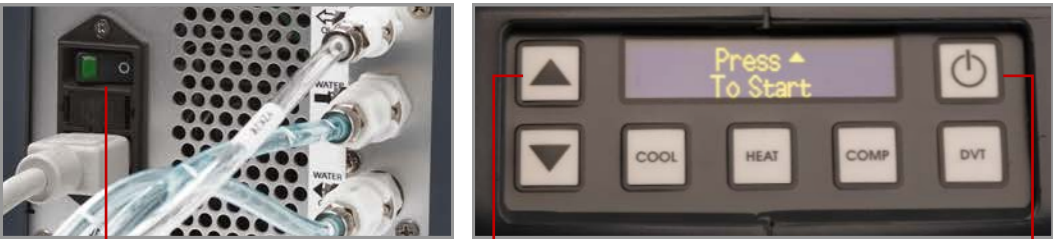


- STEP 2

Flip the switch *on* at the back of the machine.
- STEP 4

Press **▲** button to start therapy.
- STEP 3

Hold **⏻** button down until screen *turns on*.

4 STOPPING THERAPY

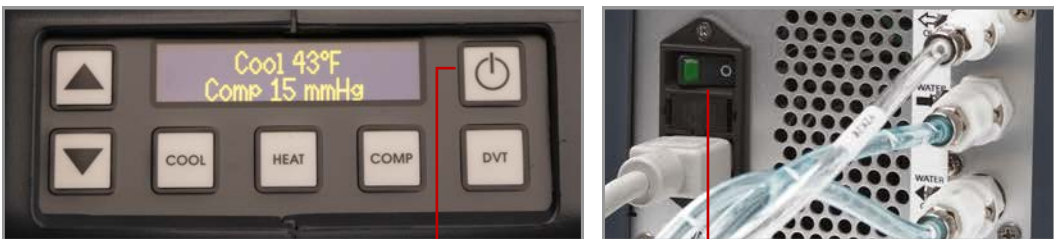

- STEP 1

Hold **⏻** button down until screen *turns off*.
- STEP 2

\*Flip the switch *off* at the back of the machine.  
\*Important: Do not switch off in the back before turning off the front of the unit

5 APPLICATION OF FOOT/ANKLE COLD COMPRESSION & DVT WRAPS


**STEP 1** Place single DVT wrap around opposite calf of surgical leg, hose going down towards foot.



**STEP 2** Connect umbilical hose to cold compression wrap.



**STEP 3** Fold and velcro long end of wrap together to form boot. Place foot in the boot and velcro over shin then velcro over top of the foot.

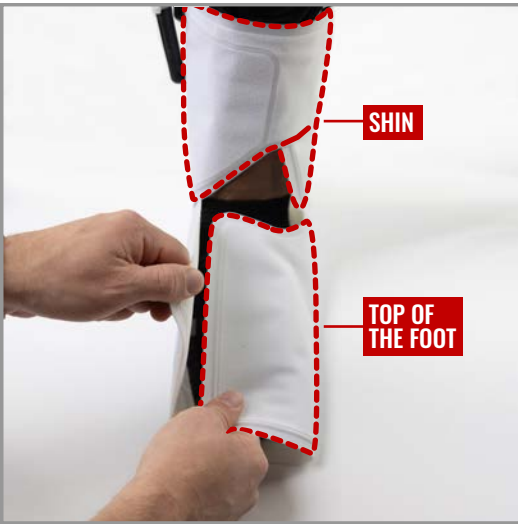

COMPLETED FOOT/ANKLE SET UP



**WATCH HOW TO SETUP THE ORTHOSAVE®**  
Scan this code on your smartphone to watch a video demonstrating steps 2-5.

6 TROUBLESHOOTING

MACHINE / SYMPTOM	POSSIBLE SOLUTION(S)
<b>ORTHOSAVE® DVT</b> The OrthoSave® is beeping, but <i>shows no error code</i> on the display screen.	<ul style="list-style-type: none"><li>a. Remove DVT unit by pressing black plastic release button on the corner of the unit.</li><li>b. Shut DVT unit off then back on to reset.</li><li>c. Check connections and make sure they are plugged into the DVT ports.</li><li>d. Check for pinched wrap air lines.</li><li>e. Make sure DVT wraps are placed on the outer side of the calf with the clear tubes facing the outer side of leg/ankle.</li></ul>
<b>ORTHOSAVE® COMP</b> The OrthoSave® is beeping and <i>shows an error code</i> on the display screen. Fault codes will appear.	<ul style="list-style-type: none"><li>a. Fault code "2060": Low liquid in reservoir or hose not completely connected to wrap.<ul style="list-style-type: none"><li>▪ Shut comp unit off then turn back on to reset. Check line connections, if not connected you will hear/feel snap.</li><li>▪ How to refill the liquid reservoir when low:<ul style="list-style-type: none"><li>1. Start the OrthoSave® Comp.</li><li>2. Remove cap from the top of the machine and begin refilling the machine with distilled water while the unit is running until the reservoir is full.</li><li>3. Once fill line is full, screw cap back to the top of the OrthoSave® Comp.</li></ul></li></ul></li><li>b. Call your Sales Service Rep for any additional questions.</li></ul>

7 RETURNING THE ORTHOSAVE®



**STEP 1** Place the UPS return label on the box, covering any previous shipping labels.

**STEP 2** Disconnect all hoses from the OrthoSave® Comp.

**STEP 3** Throw away all wraps and pads in a trash can.

**STEP 4** Place the machine, power cord, and umbilical hose in the box for return.  
(Note the separate compartments for the OrthoSave® machine, power cord and umbilical hose)

**STEP 5** Tape the box up so the machine is secure.

**STEP 6** Connect with your Sales Service Rep and have them schedule UPS to pickup the equipment.



**WATCH HOW TO RETURN THE ORTHOSAVE®**  
Scan this code with your smartphone to watch a video demonstrating on how to safely pack and return the OrthoSave®.